

A black and white photograph of a man with a prosthetic left leg running on a beach. He is wearing a dark t-shirt and shorts. A dog is running alongside him. The background shows a beach and a large, rocky cliff. The word 'ottobock.' is written in blue in the top right corner. There are two overlapping circles, one blue and one green, on the right side of the image. The text 'Quality for life' is written in white inside the green circle.

ottobock.

Quality for life

## Pure Health.

Along a beach, over a trail, on the track... there is nothing like the satisfaction and exhilaration of a good workout!

The new 3S80 Fitness Prosthesis makes it possible to get up and moving again. The only fitness system for above-knee amputees, the 3S80 can help your patients increase strength and mobility while enhancing their overall health and well-being.

Get your active patients back in action again.

# Fitness Prosthesis

with the new 3S80 Fitness Knee Joint




NEW

## 1 3S80 Fitness Knee

Able to withstand the impact of running, the knee's rotary hydraulic controls flexion and extension for dynamic response to different running speeds. The knee's rotary hydraulic can accommodate rapid flexion and extension; it won't slow the user down and will still provide a secure foot strike for powerful strides.


The 3S80's compact, lightweight design makes it ideal for either weekend jogs or competitive sprints. And patients appreciate the manual lock for security during warm-ups, stretching or other activities.

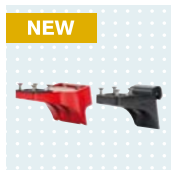

NEW

## 2 4R204, 4R206 Fitness Adapters

The 4R206 (red) allows you to quickly and easily align the foot and determine the proper height of the foot prior to cutting the pylon. Once that is established, the 4R204 (black) is used for the definitive fitting.

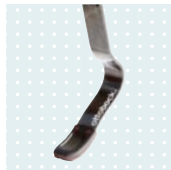

## 3 1E90 Sprinter Foot

The lightweight, durable carbon spring provides excellent performance thanks to a particularly high energy return. The 1E90 is available in six different stiffness levels to provide the precise blend of dynamic performance depending upon patient indications and desired running styles.



NEW

## 4 2Z500, 2Z501 Fitness Soles

Two new soles add even more versatility and performance to the 1E90 Sprinter. The all-terrain version (2Z500) provides just the right amount of grip for a variety of running surfaces while the spiked sole (2Z501) is designed specifically for use on a track.



## Try it! Ask for the 3S80 Fitness Prosthesis Trial Kit

Now your patients can experience the 3S80 Fitness Prosthesis. Ottobock's unique trial kit includes the 3S80, a variety of adapters and a series of 1E90 feet to help you find the perfect balance of fitness prosthetic components for your patient. Contact your Ottobock Sales Representative to schedule a trial fitting. If you know which products your patients need, simply call our Customer Service at 800 328 4058.

## Technical Data 3S80 Fitness Knee

Max. body weight	220 lbs
Weight	682 g
Max. flexion angle	135°