

ottobock.



Quality for life

Zenit CLT

Easy handling with an aesthetic design

This folding wheelchair, which has the look and driving characteristics of a rigid-frame wheelchair, stands out for its low weight and its styling. The centre of gravity, back support and caster wheel connection of the open-frame wheelchair are permanently welded. Thanks to the innovative Easy-Go folding system and the foot support that also folds up, the Zenit CLT is easy to transport. The Zenit CLT is available as an aluminium version in various frame colours or in carbon.

Benefits at a glance

- Stylish high-end folding wheelchair
- Permanently welded yet foldable
- Configurable, tailor-made wheelchair
- Open-frame construction with permanently welded centre of gravity, caster wheel attachment and back support
- Straightforward to fold and unfold with the Easy-Go folding system and the foot support that also folds up




- Permanently welded frame lends the driving characteristics and look of a rigid-frame wheelchair.



- For fine-tuning: more active or passive centre of gravity adjustment by 12 mm.



- The Zenit CLT can be ordered with a flexibly adjustable or permanently welded lower leg length.

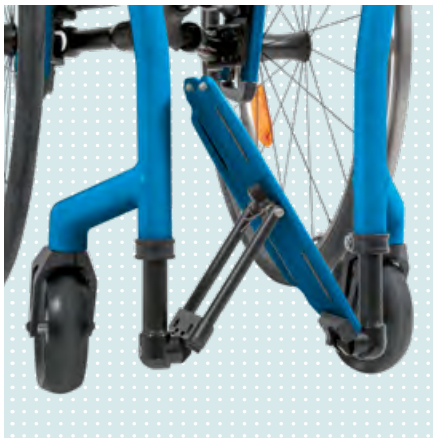

- Another special feature of the Zenit is the automatically folding foot support.



- The back support pad of the Zenit CLT is flush with the ergonomic push handles.

Technical data

Seat width	320–460 mm
Seat depth	360–500 mm
Front seat height	380–540 mm
Rear seat height	370–500 mm
Back support height	250–500 mm
Lower leg length	Adjustable: 280–510 mm Permanently welded
Anterior frame angle	75° / 85°
Overall width	490–725 mm
Overall length	720–1.040 mm
Overall height	720–1.040 mm
Max. load	100 kg (220 lbs)
Weight from approx.	8.9 kg

Powder coating



Sparkle light orange



Vinho sparkle



Candy red



Marine blue




Sky blue



Sparkle granny-smith



Ice blue

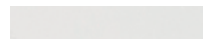

White RAL 9016



Brilliant silver




Anthracite metallic



Shimano matte



Jet black RAL 9005

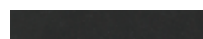

Jet black RAL 9005
matte (texture)

The additional style package includes coloured drive wheel rims and caster wheel adapters in the colours sparkle light orange, marine blue, candy red and sparkle granny-smith.

Standard model

- Permanently welded centre of gravity that can be adjusted by 12 mm
- Easy-Go folding system
- Weight 8.9 kg and up

Options

- Carbon side panels
- Carbon foot support
- Various drive wheels and handrims
- Various back support and seat upholstery
- Style package

