


ottobock.

28U25 WalkOn Reaction junior

It's so dynamic!

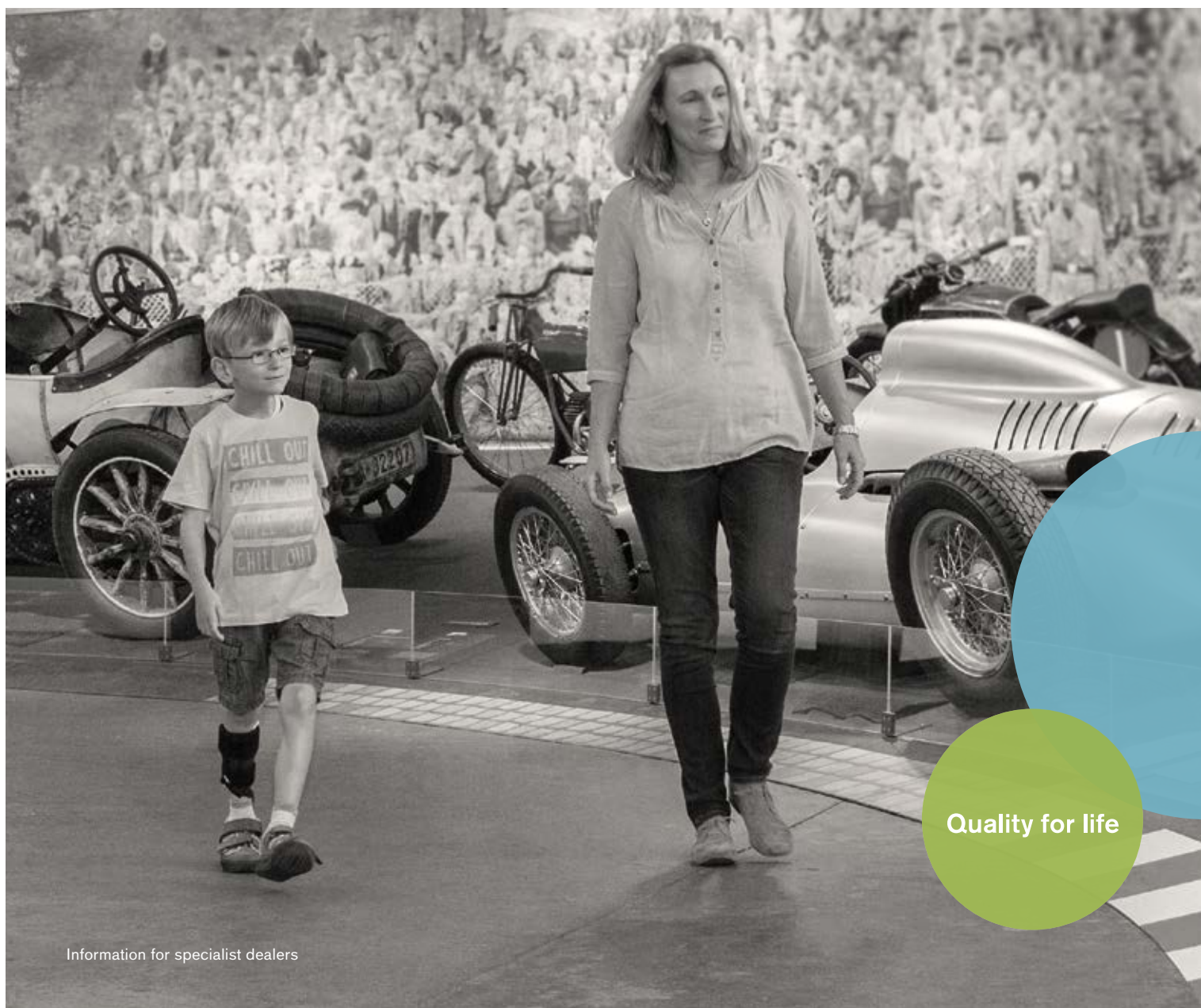

Quality for life

WalkOn Reaction junior

Dynamic ankle foot orthosis
made of carbon for children

The WalkOn Reaction junior is the new member of the WalkOn line. Based on the WalkOn Reaction plus, it features highly dynamic properties. The ankle foot orthosis uses the ground reaction force to significantly improve the user's gait pattern. It is geared towards children, who require greater support than that provided by traditional dorsiflexion assist orthoses.



Benefits at a glance

- Highly dynamic properties (use of ground reaction force)
- Carbon fibre prepreg material for high energy return
- High wearer comfort due to low weight and open heel and the anatomically designed frontal support element
- Shin pad made of climate-regulating material for improved skin hygiene and wearer comfort
- Combination with additional dynamic ankle foot orthoses (DAFOs) possible
- All fabric parts can be washed and replaced
- Sole can be trimmed to size
- Can be worn with sturdy standard shoes

Indications

- Cerebral palsy
- For dorsiflexor weakness, without or with mild to moderate spasticity
- Hemiparesis, diparesis

Effects

- Supports lifting of the foot
- Promotes dynamic rollover
- Increases the stride length by improving muscle dynamics
- Promotes exterior rotation of the foot at heel strike
- Promotes a physiological gait pattern
- Reduces the risk of falling





The orthosis design is adapted to the anatomical conditions of children and for this reason differs in its various sizes. Because it is comfortable to wear, the WalkOn Reaction junior enjoys high acceptance amongst children.

Art. no. 28U25

Order no.: Article number=side & shoe size | Example for ordering: 28U25=L24-27



Side	Shoe size	Overall height (cm)
L/R	24-27	20.5
L/R	27-30	24.5
L/R	30-33	27
L/R	33-36	29.4

A FIT KIT (28T5) is available for the WalkOn Reaction junior, as it is for the other products in the WalkOn line.

© (2016) Tabaluga Enterprises GmbH.
TABALUGA is based on an idea from
Peter Maffay, Gregor Rottschalk and
Rolf Zuckowski. Images by Helme Heine.
www.tabaluga-enterprises.de

Otto Bock HealthCare GmbH
Max-Näder-Straße 15 · 37115 Duderstadt/Germany
T +49 5527 848-1706 · F +49 5527 848-72330
export@ottobock.de · www.ottobock.com